



GoLift
Product Innovation
Award Winner!

Amico




The Amico GoLift is setting new standards in aesthetics, performance, durability, ease of use and ease of service.



GoLift Casework and Tracking for every need



GoLift tracking includes curves, turntables, XY gantries, gantry gates, and is readily enabled with constant charge and can be ceiling mounted or recessed.

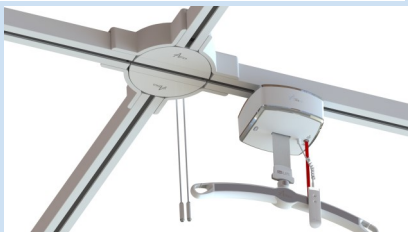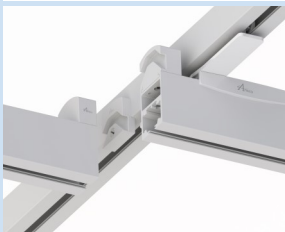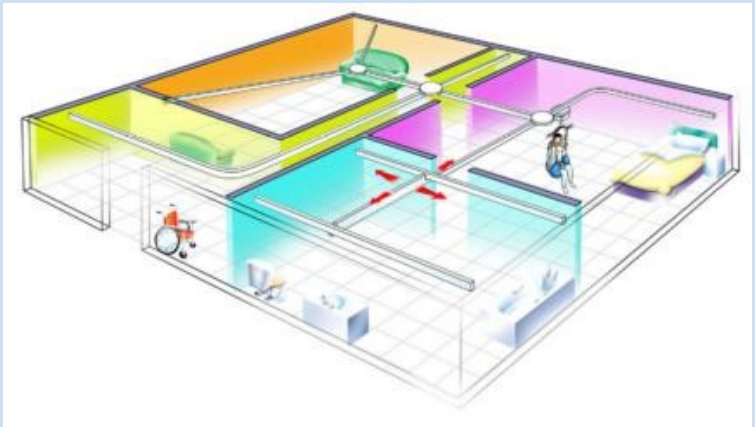

Update your lifter with **GoLift** not your track

GoLift is equipped with an easy to attach trolley which allows for instant mechanical and electrical connection of the **GoLift**. This allows effortless moving of the **GoLift** without tools when servicing or moving from track to track. This allows for a facility to save cost by having to purchase fewer lifts to cover the entire facility. **GoLift** have purpose designed trolleys for most major brands enabling upgrade and replacement without changing tracks. Facilities with tracks from different brands will be able to homogenous to one hoist and easily move about as desired.



GoLift excels at aesthetics and infection control



GoLift ceiling lift, carry bar and hand control incorporate smooth edges and rounded corners for infection control.

GoLift delivers uncompromised quality and performance while being aesthetically pleasing and the most compact ceiling lifter/hoist in its class.

GoLift incorporates an impressive set of all metal gears and state-of-the-art battery technology that allow the caregiver to safely, and effortlessly, transfer a patient weighing up to 317kgs (single cassette). The long life NiMH batteries can be easily replaced when necessary without tools and without the need for a service technician.

GoLift lifting strap is aligned along the track axis minimizing wear and tear and maximizing strap life.



GoLift
Product Innovation
Award Winner!

*A*mico

Features of *GoLift 400, 700, 1000*

- **GoLift400** (181 kg) / **GoLift700** (317 kg)
- 197mm L x 197mm W x 108mm H
- **GoLift1000** lbs (453 kg)
(using two GoLift 700 lifts in tandem)
- 394mm L x 197mm W x 108mm H
- Emergency Stop/Lowering Device
- Upper /Lower Limit Detection
- Slack Tape Sensor
- Over speed governor
- Low Battery and Dead Battery Alarms
- Soft Start and Stop
- Overload Protection
- Multiple charging options including continuous charge
- Case Flame retardant ABS
- Sound <65 dB
- Unit weight 8, 10, 20 lbs
- Protection Class: IPX0
- Capacitive touch hand control eliminating hand stresses
- Hand control attaches magnetically to the carry bar. It never dangles in the way
- Operating and battery warning lights are on the hand control maximizing convenience and visibility.
- Charger 100-240 VAC 40 Watt
- Charging time: 2-4 hours
- Lifting speed 2"/sec loaded 1" /sec
- High Capacity Ni-Mh battery
- Std: 2 x 14.4V (2Ah) Opt: (3.3 or 5 Ah)
- Straight or curved tracks
- X-Y, Gates, Turntable
- Recessed "hidden" track

Rehab Assist
Making Life Better

Tel 03 97297978 / 1300364534
sales@rehabassist.com.au
www.rehabassist.com.au
PO Box 536 Bayswater Vic 3153