

WE'RE THE COMPANY THAT BRINGS YOU BRIAN...




AND BRIAN BRINGS YOU THE PRESSURE FIT RAIL SYSTEM

Introducing, the unique and portable Waverley Glen Pressure Fit System. With each component weighing no more than 25lbs, the Pressure Fit System is the newest in temporary ceiling lift solutions. Minimal installation is required because the system works by exerting equal pressure between the floor and ceiling to stabilize the unit. Constructed of high strength aluminum with a 440 lb weight capacity the Pressure Fit is suitable for most situations. Our unique patent pending one step height adjustment allows for set up in less than two minutes and the incorporated level and pressure gauge indicator allow for easy and accurate installation.

WAVERLEY GLEN
A PRISM MEDICAL COMPANY

U.S.A. 45 Progress Parkway, Maryland Heights, MO USA 63043 | T 314.692.9135 | 1.866.891.6502
Canada 87 Sharer Road, Vaughan, ON Canada L4L 8Z3 | T 905.850.0093 | 1.800.265.0677
www.waverleyglen.com

WAVERLEY GLEN Pressure FIT RAIL SYSTEM

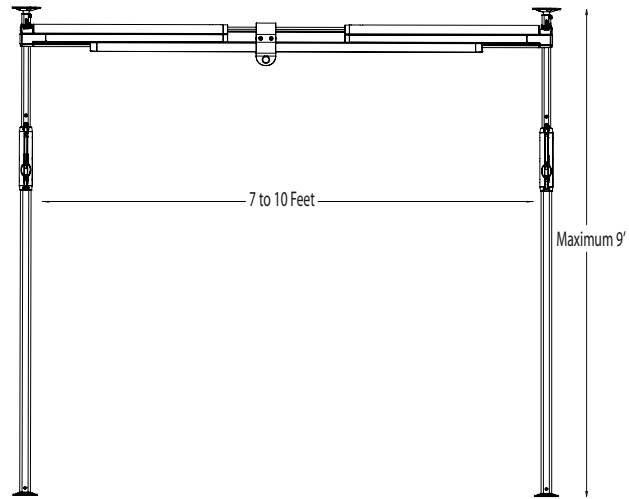

WE'RE THE COMPANY THAT BRINGS YOU GREAT THINGS IN PORTABLE PACKAGES.

PRODUCT HIGHLIGHTS

- Patent pending one step height adjustment – set up in less than two minutes
- Constructed of High strength aluminum with a 440 lb weight capacity
- Height adjustable from 7'-9'
- Slim profile uprights – approximately 1" wide
- Track adjustable in length from 7' - 10'
- Each component weighs no more than 25lbs
- Incorporated level and pressure gauge indicator for accurate installation
- Pressure between floor and ceiling will ensure that ceilings are not damaged
- Optional carry bag

TECHNICAL SPECIFICATIONS

- Width: 7' minimum, 10' maximum
- Height: Up to 9'
- Unit Weight: under 45 lbs (20 kgs)
- Maximum Load: 440 lb (200 kgs)

WAVERLEY GLEN
A PRISM MEDICAL COMPANY

Waverley Glen Systems Ltd.

U.S.A. 45 Progress Parkway, Maryland Heights, MO USA 63043 | T 314.692.9135 | 1.866.891.6502

Canada 87 Sharer Road, Vaughan, ON Canada L4L 8Z3 | T 905.850.0093 | 1.800.265.0677

www.waverleyglen.com